



PCT/GB2004/005136



INVESTOR IN PEOPLE

The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

REC'D 13 JAN 2005

WIPO

PCT

PRIORITY DOCUMENT

SUBMITTED OR TRANSMITTED IN
COMPLIANCE WITH RULE 17.1(a) OR (b)

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed

Andrew Gersy

Dated

29 December 2004





THE PATENT OFFICE
C

13 DEC 2003

NEWPORT

15AEP03 5859224-1 D02819
0.00-0328912.1 CHEQUE

The Patent Office

Cardiff Road
Newport
South Wales
NP10 8QQ

Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)

1. Your reference

MJL/C277/0

2. Patent application number

(The Patent Office will fill this part in)

0328912.1

3. Full name, address and postcode of the or of each applicant (underline all surnames)

Oxoid Limited
Kingsland Business Park, Wade Road
Basingstoke
Hampshire RG24 8PW

Patents ADP number (if you know it)

647875001

If the applicant is a corporate body, give the country/state of its incorporation

UK

4. Title of the invention

Medium for Identification of
Micro-organisms

5. Name of your agent (if you have one)

Keith W Nash & Co

"Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)

90-92 Regent Street
Cambridge CB2 1DP

Patents ADP number (if you know it)

1206001✓

6. Priority: Complete this section if you are declaring priority from one or more earlier patent applications, filed in the last 12 months.

Country

Priority application number
(if you know it)

Date of filing
(day / month / year)

7. Divisionals, etc: Complete this section only if this application is a divisional application or resulted from an entitlement dispute (see note f)

Number of earlier UK application

Date of filing
(day / month / year)

8. Is a Patents Form 7/77 (Statement of inventorship and of right to grant of a patent) required in support of this request?

Answer YES if:

- a) any applicant named in part 3 is not an inventor, or
- b) there is an inventor who is not named as an applicant, or
- c) any named applicant is a corporate body.

Yes

Otherwise answer NO (See note d)

Patents Form 1/77

9. Accompanying documents: A patent application must include a description of the invention. Not counting duplicates, please enter the number of pages of each item accompanying this form:

Continuation sheets of this form

Description 14

Claim(s) 2

Abstract 1

Drawing(s) 1 + 1 *RM*

10. If you are also filing any of the following, state how many against each item.

Priority documents -

Translations of priority documents -

Statement of inventorship and right to grant of a patent (Patents Form 7/77) -

Request for a preliminary examination and search (Patents Form 9/77) 1

Request for a substantive examination (Patents Form 10/77) -

Any other documents (please specify) -

11. I/We request the grant of a patent on the basis of this application.

Signature(s)

Keith W. Nash & Co. Date 12/12/2003
KEITH W NASH & Co

12. Name, daytime telephone number and e-mail address, if any, of person to contact in the United Kingdom

M J Lipscombe (01223) /355477

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 08459 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s). Any continuation sheet should be attached to this form.
- If you have answered YES in part 8, a Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- Part 7 should only be completed when a divisional application is being made under section 15(4), or when an application is being made under section 8(3), 12(6) or 37(4) following an entitlement dispute. By completing part 7 you are requesting that this application takes the same filing date as an earlier UK application. If you want the new application to have the same priority date(s) as the earlier UK application, you should also complete part 6 with the priority details.

C277.00/O

Title: Medium for Identification of Micro-organisms

Field of the Invention

The present invention relates to a medium for the growth and/or identification of micro-organisms, and to a method of identification of a micro-organism using the medium.

Background of the Invention

The detection and identification of micro-organisms is very important in several fields, such as diagnosis of infectious disease, tracing sources of outbreaks of disease or contamination, detecting and tracing the spread of antibiotic resistance etc.

Many methods of detecting and identifying micro-organisms have been derived in recent years, employing sophisticated techniques, such as PCR and the like. These are frequently rapid and extremely sensitive. However, there is still a need for more traditional methods involving the culturing of organisms on media (especially solid media). These traditional methods have the great advantage that they result in isolation of organisms of interest in viable form, which can then be sub-cultured into other media and/or further analysed.

It is well-known to include in such media compounds which may allow or promote the growth of certain micro-organisms whilst preventing or inhibiting the growth of other micro-organisms. Such media are said to be "selective". For example, certain antibiotics may be incorporated into media to allow the growth of micro-organisms which are sufficiently resistant to the antibiotic(s) employed.

It is also well known to include in media (which may or may not be selective) substances which allow different (possibly closely-related) micro-organisms to be distinguished. Such media may be described as "differential media".

One such differential medium is disclosed in US Patent No. 5,716,799 (Rambach). The medium disclosed in that document comprises: at least one chromogen which is a substrate for an enzyme of the particular micro-organism to be identified, together with a carbohydrate at high concentration (10-30gms/litre). The medium is such that, in the presence of the micro-organism of interest, the chromogen is hydrolysed to release a chromophore with a "derived colour" which is different to the basic colour of the chromophore in a "standard" medium. (A "standard" medium is defined in that document as being "any ordinary identification medium in which the carbohydrate has a simple function of a carbon source, at very low or even zero concentrations".)

Thus, the teaching of US 5,716,799 is that it is an essential feature of the invention described therein that the medium contains a high concentration (i.e. 10-30gms/litre) of carbohydrate. A medium in accordance with US 5,716,799 is commercially available from Becton Dickinson, USA and is known as CHROMagar™ *Candida*.

The medium disclosed in US 5,716,799 is useful in discriminating between different yeasts of the genus *Candida*. In particular, there is a clinical significance to be able to distinguish between *Candida albicans* and other species of *Candida*. *C. albicans* is a frequent human pathogen, but other *Candida* species may also be present. *Candida* species have different susceptibilities to antifungal therapy, so it is important to determine if an organism of the genus *Candida* is *C. albicans* or some other *Candia* species.

US 5,534,415 (Orenga) also discloses an identification medium, specifically a medium for the selective detection of *C. albicans*. The medium disclosed therein comprises a chromogenic or fluorogenic substrate capable of being hydrolysed by a hexosaminidase enzyme "associated" with *C. albicans*, and at least one hexosamine-containing "activator" compound (different to the substrate) present at about 1gm/litre. A medium in accordance with US 5,534,415 is commercially available (from bioMérieux, France) and is known as *Candida* ID.

Willinger *et al* (2001 J. Clin. Microbiol. 39, 3793-3795) conducted a trial to compare the performance of CHROMagar™ Candida and Candida ID, using nearly 600 clinical specimens. On CHROMagar™ Candida plates, colonies of *C. albicans* were green, whilst other colonies of other *Candida* species were pink, violet or white. On Candida ID plates, *C. albicans* colonies appeared blue, whilst other *Candida* species were pink or white.

Willinger *et al* found that, for about 50% of specimens, identification of *C. albicans* was more rapid on Candida ID plates than on CHROMagar™ Candida plates, the blue colour on Candida ID plates typically appearing after 24 hours' incubation at 35°C, whilst the green colour of *C. albicans* colonies on CHROMagar™ Candida plates often was not apparent until incubation for about 48 hours. This was expected, in that the Candida ID medium contains a hexosamine activator which, it is presumed, accelerates induction of the hexosaminidase enzyme which acts on the chromogen, whilst CHROMagar™ Candida medium does not contain any such activator. However, they also found that "Candida ID is not as good as CHROMagar™ Candida in detecting polyfungal specimens".

In addition, CHROMagar™ Candida is superior to Candida ID in terms of its ability to allow different *Candida* species to be distinguished. Candida ID medium is very effective at distinguishing between *C. albicans* and other species of *Candida* but cannot, for example, allow for efficient differentiation between *C. tropicalis* and other non-*albicans* species. CHROMagar™ Candida, in contrast, allows for the identification of a number of *Candida* species, including *C. tropicalis*, since various species give different "derived" colour changes on the medium. Candida ID medium does not give a "derived" colour chromophore which is different to that obtained with a basic medium.

The present invention is concerned with media for micro-organisms, especially identification media and, in particular, a medium which allows for discrimination between different species of *Candida*.

Summary of the Invention

In a first aspect the invention provides a medium for the detection and/or identification of a *Candida* yeast, especially *C. albicans*, the medium comprising: a chromogen; carbohydrate in the range 1-5gms/litre; and an alcohol; the medium being such that growth of the *Candida* yeast under appropriate conditions results in hydrolysis of the chromogen to generate a chromophore of a derived colour which is a different colour from that generated by hydrolysis of the chromogen in a standard medium.

For present purposes, a "derived colour" means any colour (generated by hydrolysis of the chromogen in a medium in accordance with the invention) whose dominant wavelength differs from the dominant wavelength of the chromophore generated by hydrolysis of the chromogen in a "standard medium". Again for the purposes of the present specification, a "standard medium" may be defined as a medium essentially lacking an alcohol but otherwise in accordance with the invention. In particular, the derived colour is preferably produced by hydrolysis of the chromogen as a result of the growth and/or presence of *C. albicans* in/on the medium, such that the presence of *C. albicans* may be determined by the appearance in the medium of the chromophore exhibiting derived colour.

The dominant wavelength of the chromophore may be determined by reference to daylight, as defined by the CIE, using any conventional technique for measuring the colour of an object (e.g. with a spectrophotometer).

In most instances the colour difference between the derived colour of the chromophore and the normal colour of the chromophore will be readily apparent to the naked eye of a human observer.

The medium is especially useful in the detection and/or identification of *C. albicans* and *C. tropicalis*, and in particular in allowing differentiation between *C. albicans* and other *Candida* species. The medium is preferably a solid medium.

In a second aspect the invention provides a method of detecting and/or identifying a *Candida* yeast in a sample (especially wherein the sample may contain other micro-organisms, such as yeasts and, in particular, a plurality of other yeasts of the genus *Candida*). The method involves use of the medium of the first aspect of the invention.

More specifically, the method of the invention comprises the steps of: contacting the sample with the medium of the first aspect of the invention defined above, under appropriate conditions, to allow growth of the *Candida* yeast; and detecting the presence of a chromophore having a colour indicative of the presence of the *Candida* yeast. The method of the invention is especially useful for, but not necessarily limited to, distinguishing *C. albicans* from other *Candida* species.

The inventor has surprisingly found that, contrary to the teaching of US 5,716,799, it is possible to discriminate between *Candida* species using a medium which does not contain a high concentration of carbohydrate.

Carbohydrate

Carbohydrate used in the medium of the present invention may, in principle, be any carbohydrates which are capable of being used as a carbon source by *Candida*. The carbohydrate content preferably comprises at least some sugars, especially mono- or disaccharides. In particular pentose or hexose monosaccharides are especially convenient. The most preferred sugar is glucose. The carbohydrate content may additionally or alternatively comprise more complex saccharides with a higher degree of polymerisation. A convenient source of carbohydrate is malt extract, which is readily available commercially. In one embodiment the medium comprises both glucose and malt extract in an amount such that the total concentration of carbohydrate (glucose plus malt extract and those carbohydrates (if any) present in the peptones) in the medium is in the range 1-5gms/litre.

The carbohydrate is preferably incorporated into the medium at a concentration in the range 2-4gms/litre.

Alcohol

The alcohol is preferably present in the medium at a concentration which does not exert a significant inhibitory effect on the growth of *Candida* yeasts. A suitable concentration for most alcohols is in the range 1-10mls/litre, preferably 2-8, more preferably 5-7mls/litre.

The alcohol may be an aromatic, alicyclic, or aliphatic (straight chain or branched) compound, although aliphatic compounds are generally preferred.

More especially the alcohol preferably comprises between 1 and 5 carbon atoms, preferably 2-4 carbon atoms. Ethanol (99.7-100% v/v) is presently the preferred alcohol, although industrial methylated spirits ("IMS"), comprising a mixture of about 85% ethanol and about 10% methanol, is also suitable and is commercially available at a cost lower than pure ethanol.

Chromogen

The medium of the invention comprises at least one chromogen, although two or more chromogens may be used if desired. The chromogen is a compound which, in the presence of species specific enzymes and in suitable conditions, is hydrolysed to release a chromophore which is discernibly different to the chromophore generated by growth of the same yeast in a standard medium containing the same chromogen.

The chromogen preferably comprises a hexosaminide, most preferably a glucosaminide, and typically conforms to the general formula X-GluNAc (wherein Glu is glucose, NAc is an N-acetyl group and X is a chromophore). Many suitable chromogens are known to those skilled in the art. A preferred chromogen is 5-bromo-4-chloro-3-indolyl N-acetyl β -D-glucosaminide. Other chromogens which may additionally or alternatively be used include 5-bromo-6-chloro-3-indolyl N-acetyl β -D-glucosaminide, X-Gal NAC (wherein Gal is galactose, NAC is an N-acetyl group and X is a chromophore), 5-bromo-6-chloro-3-indolyl phosphate *p* toluidine salt, 5-bromo-4-chloro-3-indolyl phosphate *p* toluidine salt and 6-chloro-3-indoxylphosphate.

Other Components

The medium will generally comprise additional components which are present in a typical conventional yeast medium, such as agar, yeast extract, sodium chloride and one or more antibiotics (e.g. chloramphenicol) to suppress bacterial growth.

In a preferred embodiment, the medium may also contain one or more of the following: malic acid (at about 0.1gms/litre); peptones (at about 2.5gms/litre); and a buffer (e.g. KH_2PO_4 , at about 2.0gms/litre).

The method of the invention involves contacting a sample to be tested with a medium in accordance with the first aspect of the invention. Typically the sample will comprise a suspension of yeast cells and other material and may be, for example, a clinical sample obtained from a patient. The sample will normally be spread over the surface of a solid medium, contained within a Petri dish or other suitable container. The medium to which the sample has been applied must then be incubated at an appropriate temperature and for sufficient length of time to allow growth of any yeast cells in the sample to a sufficient extent to cause hydrolysis of the chromogen. Typically the temperature of incubation is in the range 10-37°C, more preferably 20-35°C, and most preferably 25-30°C. The length of incubation depends, at least in part, on the incubation temperature, but for an incubation temperature in the range 25-30°C an incubation of 24-48 hours duration is usually adequate. If a result is not apparent after 24 hours, the plates can simply be returned to the incubator and re-inspected at say, 36 or 48 hours.

Typically, for a medium in accordance with preferred embodiments of the invention, the colour of the chromophore will be green or green/blue in the presence of *C. albicans*, purple or blue in the presence of *C. tropicalis* and pink in the presence of *C. krusei*, where the chromogen is 5-bromo-4-chloro-3-indoxyl N-acetyl β -D-glucosaminide. Green and green/blue constitute derived colours: the normal colour of the chromophore in a standard medium is purple or blue. These colour differences are readily apparent to the naked eye.

If desired, control cultures of known *C. albicans* and/or known *C. tropicalis* strains, for example, can also be set up to give a comparison with the test culture.

For the avoidance of doubt, it is expressly stated that any feature of the invention described as “preferred”, “advantageous”, “desirable” or the like may be employed in any embodiment of the invention in isolation or in combination with any other feature or features so described, insofar as the context permits.

The invention will now be further described by way of illustrative example.

Examples

Example 1

A basic yeast medium formulation was prepared, according to the Table below:

Table 1 – Basic yeast medium

Ingredient	gm/l
Agar	15
Yeast extract	2
Peptone	2.5
Malic acid	0.1
Glucose	0.5
Malt extract	1.5
Sodium chloride	2.5
5-bromo-4-chloro-3-indoxyl N-acetyl β -D glucosaminide	0.1
5-bromo-6-chloro-3-indoxyl phosphate <i>p</i> toluidine salt	0.08
Chloramphenicol	0.1
(Total	24.38 gm/l)

Various test media (in accordance with the invention) or comparative media (outside the scope of the invention) were prepared by supplementing the basic formulation as follows:

Table 2

Example	Basic formulation supplemented with:	Notes
A	6ml IMS	1
B	6ml ethanol	2
C	6ml butanol	
D	CHROMagar™ Candida	3 6
E	4ml methanol	4
F	6ml propan-2-ol	
G	2gm/l sodium succinate	6
H	2gm/l sodium citrate	6
I	6ml/l ethyl acetate	5 6
J	6ml/l ethanol	5
K	6ml/l H ₂ O	5 6

Notes

1. IMS is approximately 85% ethanol, 10% methanol.
2. Ethanol was at least 99.7% ethanol v/v.
3. 47.7 gm/l CHROMagar™ Candida, for comparison.
4. (6ml/l methanol was found to inhibit growth of many *Candida* strains).
5. Media I, J and K contained carbohydrate at 4gm/l, twice the concentration of carbohydrate in media A-H. This was obtained simply by doubling the amount of glucose and malt extract to 1gm/l and 3gm/l respectively.
6. Comparative example, outside the scope of the invention.

The final pH of formulations A-C and E-H was adjusted to 6.0 with NaOH or HCl as appropriate. The commercially available CHROMagar™ Candida medium (formulation 'D') was found to have a pH of 6.03. The pH of formulations I-K was found to be about 5.8 and was not adjusted.

The various media formulations were tested for their ability to give rise to a derived colour change with a number of *Candida* yeast strains, which were either available within the Oxoid yeast culture collection (OYC) or were clinical isolates. The table below shows the yeast strains tested (in this or subsequent example):

Table 3 – Yeast Strains

Strain	NCPF equivalent	ATCC equivalent	Notes
OYC247 (<i>C. albicans</i>)	3153		
OYC205 (<i>C. albicans</i>)	3179	10231	NCYC1363
OYC263 (<i>C. albicans</i>)		2091	NCYC854
OYC238 (<i>C. albicans</i>)		18804	NCYC597
OYC249 (<i>C. glabrata</i>)	3240		
OYC225 (<i>C. kefyr</i>)		8555	
OYC223 (<i>C. tropicalis</i>)		750	
NCYC4 (<i>C. tropicalis</i>)			NCYC4
OYC222 (<i>C. krusei</i>)		6258	

NCPF – National Collection of Pathogenic Fungi, PHLS Mycological Reference Laboratory, 61 Colindale Avenue, London NW9 5HT, United Kingdom.

ATCC – American Type Culture Collection, Manassas, Va 20108, USA.

NCYC – National Collection of Yeast Cultures, Institute of Food Research, Norwich Research Park, Colney, Norwich NR4 7UA, United Kingdom.

All strains without an NCPF or ATCC equivalent are available on request from Oxoid Limited.

The various strains were grown on Sabouraud Dextrose agar plates. Growth from the plates was then used to prepare suspensions in saline solution, each suspension being adjusted to a turbidity of 1 McFarland Unit. Pure cultures of these suspensions were then streaked onto the various agar formulations and incubated at 30°C for up to 48 hours. The results are shown below:

Table 3A

Formulation	OYC247 (<i>C. albicans</i>)	OYC263 (<i>C. albicans</i>)	NCYC4 (<i>C. tropicalis</i>)
A	Green	Dark green	Purple-blue
B	Green	Dark green	Purple
C	Green-blue	Blue-green	Blue-green
D	Green-yellow	Green-grey	Purple-blue
E	DarkBlue/Blue-green	Blue-violet/Green	Purple
F	Dark Blue/Blue-green	Violet/Blue-green	Purple
G	Dark blue/Blue-green	Blue-green	Pink-purple
H	Dark blue	Violet/Blue-green	Pink-purple
I	Blue-green	Blue-violet/Blue-green	Purple/Blue-grey
J	Green	Green	Purple
K	Blue/Blue-light green	Blue-light green	Purple/Blue-grey

Formulation	OYC223 (<i>C. tropicalis</i>)	OYC222 (<i>C. krusei</i>)	OYC 249 (<i>C. glabrata</i>)	OYC 225 (<i>C. Kefyr</i>)
A	Blue-purple-grey	Pink-purple	Not done	Not done
B	Blue-purple-grey	Pink-purple	" "	" "
C	Blue-green	Pink-purple	" "	" "
D	Blue-purple	Pink-Lt purple	" "	" "
E	Dark blue/Dark purple	Pink-purple	" "	" "
F	Purple-blue-grey	Pink-purple	" "	" "
G	Purple-Lt grey	Pink-purple	" "	" "
H	Purple-Lt grey	Pink-purple	" "	" "
I	Blue-violet	Pink	" "	" "
J	Blue-grey	Pink	Light pink/beige	Dark pink
K	Blue-purple	Pink	Beige	Pink

From the results it is clear that only the agar plates comprising the combination of alcohols and low carbohydrate concentration were able to give a derived colour change for *C. albicans*. Other carbon sources, such as succinate (formulation G), citrate (formulation H) or ethyl acetate (formulation I) were ineffective.

The most effective formulation in accordance with the invention, in terms of giving a clear derived colour change for *C. albicans*, was formulation J, comprising 6mls/l ethanol and 4gms/l carbohydrate. Accordingly, ethanol (and, to a lesser extent, IMS) is the preferred alcohol for use in media and methods in accordance with the present invention. Formulations A and B also gave acceptable results.

Example 2

A further experiment was performed to compare a medium in accordance with the present invention (essentially corresponding to formulation J as described in Example 1), with the commercially available media Candida ID and CHROMagarTM Candida.

Pure cultures of various *C. albicans* strains were grown on Oxoid Sabouraud Dextrose agar plates for 48 hours, and these used to prepare cell suspensions in saline solution having a McFarland turbidity equivalent to standard 1. The suspensions were then streaked onto test plates using the diminishing streak technique, and incubated at 30°C for 24 hours.

After 24 hours, colour was generally apparent only on the most heavily inoculated parts of the plates (first and second zones of inoculation). Single colonies, on all plates, generally lacked colour. The results are shown below.

Table 4

<i>C. albicans</i> Strain	Candida ID (BioMerieux Lot 776691301)	CHROMagar Candida(BD) Lot 402-c	Invention Medium
PAD 32121	Blue	Green	Green
PAD 32365	Blue	Green	Green
OYC 238*	White	White	Green
PAD 32605	Pale blue	Green	Green
PAD 31989	Blue	Green	Green
PAD 31986	Blue	Green	Green
PAD 32755	Blue	White	Green
PAD 32770	Blue	Green	Green
PAD 32651	Colourless (ppt)	Colourless (ppt)	Colourless (ppt)
PAD 32433	Pale blue	Green	Green
PAD 32617	Pale Blue	Green	Green
OYC 205*	Blue	Green	Green
PAD 32777	Blue	Green	Green
PAD 32611	Blue	Green	Green
PAD 222973	Blue	Green	Green
PAD 31987	Blue	Green	Green
PAD U32184	Blue	Green	Green
OYC 263*	Blue	White	Green
PAD 31984	Blue	Green	Green
PAD 32289	Blue	Green	Green
Total positive (24h)	18/20	16/20	19/20

PAD = Pure clinical isolates from St. Mary's Hospital, Paddington, London.

OYC = Oxoid Culture Collection

Candida ID – *C. albicans* should give a blue colour (i.e. no derived colour change)

* see Table 3

All the strains used in this experiment are available from Oxoid upon request.

The results show that, whilst all three media performed well, after 24 hours' incubation the medium in accordance with the invention returned more positive results than the CHROMagar™ Candida medium. The medium of the invention was also comparable, or possibly even slightly superior, to the Candida ID medium in terms of positive results obtained after 24 hours. This was extremely surprising and unexpected (in relation to Willinger *et al*, who found 335/384 *C. albicans* grown on Candida ID compared to 167/370 grown on CHROMagar™ Candida that were correctly identified at 24 hours), given that the Candida ID medium comprises a hexosaminide activator to accelerate induction of the relevant hexosaminidase.

After 48 hours incubation, all media gave positive results with all 20 strains tested, and most individual colonies were coloured. However, growth of PAD 32651 on the Candida ID medium was poor, and many individual colonies were colourless or only faintly coloured.

Example 3

A suspended mixture of different species of *Candida* was used to inoculate a solid medium in accordance with the invention (essentially corresponding to formulation J as described in Example 1), and incubated at 30°C for 48 hours. Typical results are shown in Figure 1, which is a photograph of a plate showing the growth of the *inoculum* and separate colonies. Even in monochrome, the difference between colonies of different species is apparent.

Claims


1. A medium for the detection and/or identification of a *Candida* yeast, the medium comprising: a chromogen; carbohydrate in the range 1-5gms/litre; and an alcohol; the medium being such that growth of the *Candida* yeast under appropriate conditions results in hydrolysis of the chromogen to generate a chromophore of a derived colour which is a different colour from that generated by hydrolysis of the chromogen in a standard medium.
2. A medium according to claim 1, wherein the chromogen is hydrolysed in the presence of *C. albicans* to give a chromophore with a derived colour.
3. A medium according to claim 1 or 2, comprising comprising carbohydrate in the range 2-4gms/litre.
4. A medium according to any one of claims 1, 2 or 3, comprising glucose.
5. A medium according to any one of the preceding claims comprising malt extract.
6. A medium according to any one of the preceding claims, comprising an alcohol in the range 1-10mls/l.
7. A medium according to claim 6, comprising an alcohol in the range 2-8mls/l.
8. A medium according to claim 7, comprising an alcohol in the range 5-7/l.
9. A medium according to any one of the preceding claims, comprising ethanol.
10. A medium according to any one of the preceding claims comprising 5-bromo-4-chloro-3-indolyl N-acetyl β -D-glucosaminide or 5-bromo-6-chloro-3-indolyl phosphate *p* toluidine salt or 5-bromo-6-chloro-3-indolyl N-acetyl β -D

glucosaminide or X-Gal Nac (wherein Gal is galactose, Nac is an N-acetyl group and X is a chromophore) or 5-bromo-4-chloro-3-indolyl phosphate *p* toluidine salt or 6-chloro-3-indoxyl-phosphate.

11. A medium according to any one of the preceding claims comprising one or more of the following: malic acid; peptones; and KH_2PO_4 .
12. A method of detecting and/or a *Candida* yeast in a sample, the method comprising the steps of: contacting the sample with a medium in accordance with any one of the preceding claims, under appropriate conditions, to allow growth of the *Candida* yeast; and detecting the presence of a chromophore having a colour indicative of the presence of the *Candida* yeast.
13. A method of detecting and/or identifying *C. albicans* in accordance with claim 12.
14. A method according to claim 12 or 13, wherein the medium is incubated at a temperature in the range 30-37°C for at least 24 hours.
15. A method according to claim 14, wherein the medium is incubated at a temperature in the range 30-35°C for at least 24 hours.
16. A medium substantially as hereinbefore described and with reference to the examples.
17. A method substantially as hereinbefore described and with reference to the examples.

ABSTRACT**Title: Medium for Identification of Micro-organisms**

Disclosed is a medium for the detection and/or identification of a *Candida* yeast, the medium comprising: a chromogen; carbohydrate in the range 1-5gms/litre; and an alcohol; the medium being such that growth of the *Candida* yeast under appropriate conditions results in hydrolysis of the chromogen to generate a chromophore of a derived colour which is a different colour from that generated by hydrolysis of the chromogen in a standard medium.

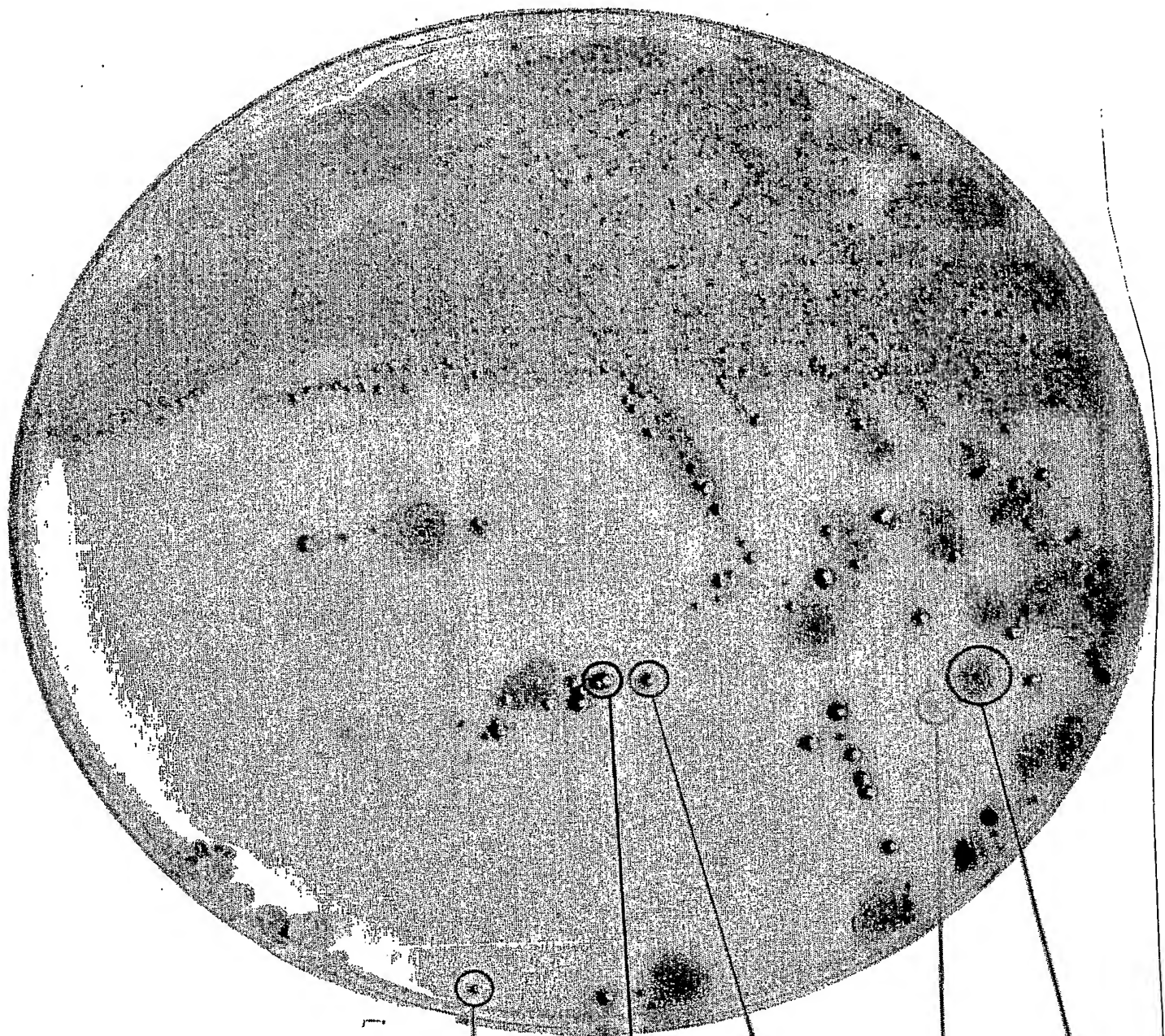

C. krusei
(fuzzy brown-pink)

C. glabrata
(beige-yellow)

C. albicans
(green)

C. tropicalis
(dark blue)

C. parapsilosis
(brown)

Fig. 1

PCT/GB2004/005136

